



RAMPART IC L148
TABLE-MOUNTED SHIELD

TABLE-MOUNTED
RADIATION PROTECTION

SHIELDING AND CURTAIN COMBINATION


The L148 Table-Mounted Shield is a system that includes abdominal shields and soft shielding curtains that provide lead-equivalencies of 1mm and 0.5mm, respectively. The L148 Multi-Purpose Shield, available separately, is an auxiliary shield with 0.5mm lead equivalency designed for versatile positioning around the patient to offer more attenuation as needed to protect both physicians and team members.

HOW IT WORKS

Working together, the Table-Mounted and Multi-Purpose Shields provide overlapping shielding designed to ensure doctors and staff are protected in a variety of different procedures.

KEY BENEFITS

- Easy Slide-On Installation to Table Rail
- Compatible With Most Cath Lab Tables
- Engineered for Durability
- 7-Part Design Provides More Protection
- 1mm Lead Equivalent Abdominal Shields
- 0.5mm Lead Equivalent Multi-Purpose Shield
- 0.5mm Lead Equivalent Soft Shielding Curtains
- Easily Replace Soft Under-Table Shielding

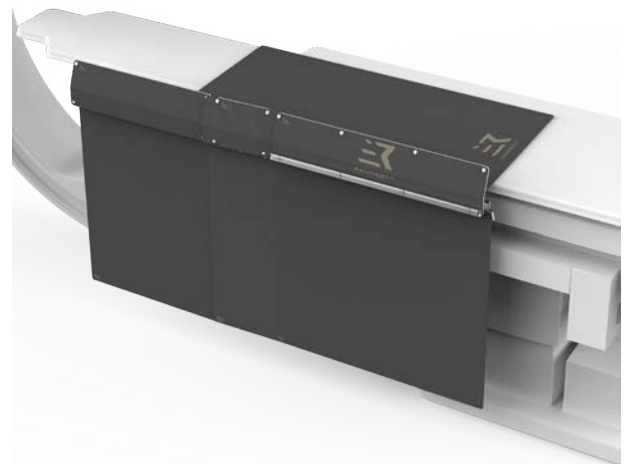

Abdominal Shields

L148 Multi-Purpose Shield (MPS)

Panel Assemblies

Adaptable Curtain Bar

MORE SHIELDING. MORE PROTECTION.

The L148 Table-Mounted Shield is made up of seven shielding components that offer more combined protective coverage than other products currently available on the market. Four curtains and three abdominal shields overlap to provide comprehensive attenuation. The L148 Multi-Purpose Shield can be positioned for additional protection as needed.

EASY ASSEMBLY AND DURABILITY

The L148 Table-Mounted Shield mounts to the procedure table by simply sliding the curtain bars onto the accessory rail. The curtain panels come pre-attached to the curtain bars, which bear the system's weight to prevent sagging of the curtains at the extreme ends. The abdominal shields easily clip onto the curtain bars.

HAZARD FREE DISPOSAL

Because the soft shield curtains come into contact with both the table and patient, over time natural wear and tear will require the curtains to be disposed of and replaced. However, no special precautions need to be taken with disposal, because our soft shielding contains no lead.



LETS TALK:

To learn more about Rampart, including a complete spec list and pricing options, please visit us at [RAMPARTic.com](https://www.RAMPARTic.com).

For ordering and pricing information, please email info@RAMPARTic.com.